

- Single Block Capacity
- Ambient +5°C to 110°C Control
- Simple to Operate
- Excellent Temperature Stability
- Wide Range of Blocks Available
- Over-Temperature Safety Cutout

### Simplicity With Stability

Ratek analogue block heaters offer excellent temperature stability combined with the ease of use of a simple analogue dial control. The DBH10 will accommodate any standard block from our extensive range, including popular 20 x 1.5mL and 12 x 15mL centrifuge tube blocks. Custom blocks in either 50mm, 80mm and 100mm depths as well as double-sided blocks are also available to cater for virtually any vessel.



### SPECIFICATIONS

DBH10	
Block capacity	1 block
Temperature control type	Solid-state analogue
Temperature control range	Ambient + 5°C to 110°C (min 28°C)
Safety features	Over temperature cut-out
Case construction	Powder-coated steel
Power input	240 VAC / 165 Watts
Overall dimensions	W165 x D175 x H90mm
Nett weight	1.6 kg (without block)

The compact design of the DBH10 makes it the ideal companion for every workstation, plus it can be extended with a nitrogen evaporation system for sample concentration, with gas manifolds available for all standard block configurations.

A custom designed solid state temperature controller is used to ensure accurate and stable temperature control, whilst the heating surface is precisely machined to ensure that excellent heat transfer and temperature uniformity is achieved. When used with Ratek machined blocks, heating is even across the block.

Once set, the DBH10 can be left unattended for extended periods of operation, with large bright heating and power indicators easily visible from across the lab.

All aluminium parts are anodised and the case powder coated for added durability and resistance to common spills.

The DBH10 is designed and made in Australia, and is backed by Ratek's unparalleled 3 year parts warranty.


### RELATED PRODUCTS

Product Code	Description
SBxx, EBxx, FBxx, CBxx	Drop-in aluminium blocks to suit virtually any vessel
EM4/6/9/12/20, EMS1	Evaporation manifolds and mounting systems to suit single block heaters

### YOU MAY ALSO LIKE

#### DBH30 - 3 Block Heater

The DBH30 is a larger analogue block heater with 3-block capacity and is ideal for working with MB1/2/3 96 well blocks.

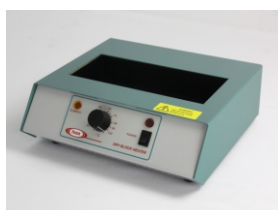

AUSTRALIAN  
MADE & SUPPORTED

3 YEAR  
PARTS WARRANTY



Quality  
ISO 9001  
SAI GLOBAL