


CAN BE
CUSTOMISED
TO SUIT YOUR
APPLICATION

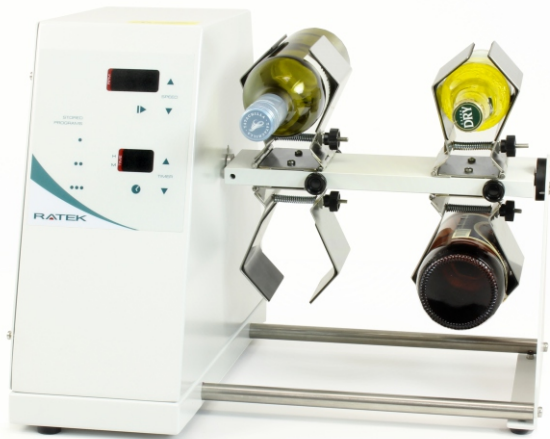

- Extremely Versatile
- Custom Made Carriers Available
- 6-57 RPM Speed Range
- End to End & Rolling Mixing Actions
- Digital Control & Large LED Displays
- Super Quiet Operation

The All-Round Rotary Mixer

The RM4 Bottle Mixer is something of an all rounder, with near endless potential for different applications and perfect for mixing wine, beer and beverage bottles.

The RM4 is ideal for the food, wine and beverage industries and delivers reproducible mixing controlled via the large LED digital control.



SPECIFICATIONS

RM4	
Mixing action	• End over end, Rolling (custom)
Capacity	4kgs evenly loaded
Speed range	6-57 RPM
Speed control type	Precision Digital + Tacho
Timer	0-99 mins or 0-99 hrs
LED speed and time display	Yes
Program presets	3 x user configurable presets
Motor and drive type	Long life super-quiet brushless AC with replaceable belt drive
Safety features	Motion interrupt cut out
Power input	240 VAC/130W
Overall dimensions	W455 x D400 x H335mm • Can be lengthened as required
Nett weight	10 kg

Thorough mixing is achieved with the end over end mixing action.

The rotating bar is mounted with precision ball bearings and is driven by a custom microprocessor controlled drive system, the large digital display can show set speed or actual speed. The in-built electronic timer can be set for 0-99 minutes or 0-99 hours, memory buttons allow for up to 3 speed / time combinations to be stored. The RM4 can be supplied in a variety of configurations to suit the application, including an inter-locked safety cover as well as a variety of clamping fixtures.

The RM4 features precision digital speed control which is complemented by the Ratek silent-belt drive system and long life brushless motor, as well as large highly visible digital LED timer and tachometer displays and audible alarm.

Built in Melbourne, supported locally and supplied with our unparalleled warranty and backup service, the RM4 is the perfect fit for a wide range of applications, and makes light work of specialised bottle mixing jobs.

ALSO AVAILABLE

Safety Enclosure

A fully enclosed guard enclosure with motion cut-out when the door is opened. Features transparent door.



RELATED PRODUCTS

Product Code	Description
UMC1/2/3/4	Universal clamps to suit RM4

AUSTRALIAN
MADE & SUPPORTED

3 YEAR
PARTS WARRANTY



Quality
ISO 9001
SAI GLOBAL